

# SIEMENS

## MODULARIS Uro Plus

**SP**

### Installation-/ start up note

Software default values

LITHOSTAR MODULARIS

© Siemens AG 1999

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

English

Doc. Gen. Date: 11.99

Print No.: SPL1-130.038.01.03.02

Replaces: SPL1-130.038.01.02.02

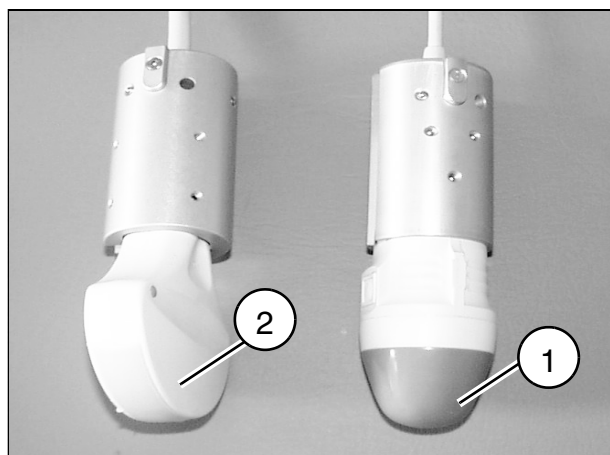
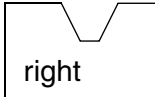


Fig. 1

	Default value
<b>Adjust c-arm turn back</b>	1200 ms
<b>install / remove option ECG-trigger</b>	not installed
<b>Pulserate / pulse rate</b>	variable1/1.5/2 Hz
<b>Ultrasound localisation</b>	<u>Sector Applicator (1/Fig. 1):</u> threshold value = 100 mm scale low = 7.7 % scale high = 7.7 %  <u>CONVEX Applcator (2/Fig. 1):</u> threshold value = 100 mm scale low = 7.7 % scale high = 7.7 %
<b>change orientation of treatment table</b>	X 

## Changes to previous version

scale high changed from 10 % to 7.7 %.

TDPS 24 / Wendt  
TDSD 24 / Groß